



Non-Invasive
Ventilation (N.I.V.)

Continuous Positive
Airway Pressure (CPAP)
and Bi-Level Ventilation.

Bi-Level Ventilation and
other NIV modes
requiring one Anti-
Asphyxia Valve

HANS RUDOLPH, inc.
8325 Cole Parkway Shawnee, KS
66227 USA
PH: (913) 422-7788
Canada: (800) 456-6695
FAX: (913) 422-3337
hri@rudolphkc.com
International:
hri@plmhealthcare.com
ISO 9001 - ISO 13485
www.rudolphkc.com

V2™ Masks Lines

There always is one V2™ Mask Version and model to effectively meet both clinical and patients needs at the Hospital, over Transportation, in Emergencies Procedures or at Home, whatever ventilator!



MULTI-PATIENT-USE / REUSABLE V2™ Line:

Steam-Autoclavable, Pasteurisable,
dishwasher safe V2™ Full Face Mask series:
Hospital and Home Care applications.

One year product servicing warranty!



SINGLE-PATIENT-USE / QUEST™ Line:

Complete QUEST™ Mask assembly can be
regularly cleaned: few mask components
easy to disassemble and to reassemble!
Ideal for Home Care applications!

Six Months product servicing warranty!



DISPOSABLE V2™ line:

Short-Term Use Disposable Full Face V2™
Mask assembly: Seven Days product
servicing warranty.

For Hospital, Emergencies and medical
transportation applications!

**A Mask Family to effectively serve Mechanical Non-Invasive
Ventilation strategies, from Hospital to Home!**

HANS RUDOLPH, inc.

V2™ Masks Lines



Three mask components only!



Reusable V2™ mask piece is made of High Grade Medical Silicone.

Quest™ Single Patient Use Mask piece and Disposable V2™ Mask piece are made of high quality soft TPE



Red and Grey Reusable Headgear for V2™ Reusable Mask series, **only Grey** Reusable Headgear for QUEST™ Mask series, **Dark Blue** Headgear for Disposable V2™ Mask series: all Headgears have 5 different adjustments and four clips with pulling tabs to allow easy and fast both mask installation and release!

MASK PIECE

- ✓ **Transparent, Latex free Lightweight**, soft material single piece made mask cone;
- ✓ **Anatomically shaped mask** with integrated inner flange to comfortably contour all faces with **Limited dead space**.
- ✓ **Chin cup** for easy mask positioning and secure wearing; **No forehead stabilizer is needed to keep mask stable and comfortable**;
- ✓ **5 different sizes** to fit more patients!

HEADGEAR

- ✓ **3 different headgears sizes**, each with 5 different adjustments;
- ✓ **Quick Easy Release Pulling tabs** on each 4 Mask Clips for fast mask removal!
- ✓ Mask assemblies deliver with headgear mounted on mask.

SWIVELS

- ✓ **Three different types of Color coded swivels** (360° + 360° rotations):

BLUE:

Integrated Silent CO2 Venting and Anti-Asphyxia Valve features. Quick-release tubing feature.

YELLOW:

Integrated Anti-Asphyxia Valve feature. Quick release tubing feature.

GREEN:

For use with ICU ventilators: these swivels integrate NO CO2 venting, NO Anti-Asphyxia-Valve. Swivels available either in 22 mm Outside or 22 mm Inside Diameters.



CO2 venting feature from single limb patient circuits.

Anti-Asphyxia-Valve:

to allow room air intake from Mask if ventilator stops for any reasons.

Quick-release Tubing Feature:

allows patients easy tubing disconnection while Mask stays on face.

HANS RUDOLPH, inc.

V2™ Masks Lines

A Mask Family to effectively serve Mechanical Non-Invasive Ventilation strategies:

From Hospital to Home!

V2™ Mask Strategies:

NIV Application Types

V2™ NIV decision based on Mechanical ventilation support length	Hospital ICU V2™ mask type	Hospital Sub-acute ICU V2™ Mask Type	Transport/emergencies V2™Mask type	Home Care V2™mask type
Short Term NIV (< or = 7 days)	Disposable V2™	Disposable V2™	Disposable V2™	Disposable V2™ & Quest™
Medium term NIV (=7days< 6 Months)	Reusable V2™	Reusable V2™		Quest™
Long-term NIV (= > 6 months)	Reusable V2™	Reusable V2™		Quest™ & Reusable V2™

Ventilators and Patient circuits Types

Type of NIV Ventilator Type of Patient Circuit	Type of RV2™ Long term NIV	Type of DV2™ Short term NIV	Type of Quest™ Medium term NIV
ICU ventilator with Expiratory valve, dual limb patient circuit:	7500 V2™	6500 V2™	6850 Quest™
Sub-Acute-ICU ventilator with Expiratory valve, single limb patient Circuit with no integrated CO2 venting valve as part of the patient circuit:	7700 V2™	6700 V2™	6870 Quest™
Home ventilator & CPAP with NO Expiratory valve, single limb patient Circuit with no integrated CO2 venting in tubing:	7600 V2™	6600 V2™	6860 Quest™

NOTE: Ventilators affording Dual limb patient circuit connection provide built-in expiratory valve

One V2™ or QUEST™ Mask Swivel to mate all circuits and to comply with all ventilators:

Don't change the mask, just replace the Swivel!

HANS RUDOLPH, inc.

V2™ Masks Lines

Mostly used part numbers - contact us for specific requirements;
ID stands for Internal Diameter (F) and OD stands for Outside Diameter (M)

REUSABLE V2™ MASK Series – One complete Mask assembly

	Large	Medium	Small	Extra-Small	Petite
7600 V2™ 22 mm OD	113485	113486	113487	113488	113489
7600 V2™ 22 mm ID	113490	113491	113492	113493	113494
7700 V2™ 22 mm ID	113500	113501	113502	113503	113504
7500 V2™ 22 mm ID	113510	113511	113512	113513	113514
7500 V2™ 22 mm OD	113505	113506	113507	113508	113509

DISPOSABLE V2™ MASK Series (one case with 18 complete mask assemblies)

	Large	Medium	Small	Extra-Small	Petite
6600 V2™ 22 mm OD	114601	114602	114603	114604	114605
6700 V2™ 22 mm ID	114701	114702	114703	114704	114705
6500 V2™ 22 mm ID	114501	114502	114503	114504	114505
6500 V2™ 22 mm OD	114506	114507	114508	114509	114510

SINGLE PATIENT USE V2™ QUEST™ MASK Series – One Complete Mask assembly

	Large	Medium	Small	Extra-Small	Petite
6860 QUEST™ 22 mm OD	113715	113716	113717	113718	113719
6870 QUEST™ 22 mm ID	113720	113721	113722	113723	113724
6850 QUEST™ 22 mm ID	113725	113726	113727	113728	113729
6850 QUEST™ 22 mm OD	113730	113731	113732	113733	113734

Key accessories V2™ Mask Lines

METRO Seal (Pack of 25):
Part n°: 669363



Mask Sizing Gauge:
Part n°: 691143



MAKERS OF RESPIRATORY VALVES SINCE 1938

- ✓ NIV to Home
- ✓ Respiratory interfaces,
- ✓ Pulmonary Function Testing,
- ✓ Stress/Exercise Testing,
- ✓ Respiratory Care,
- ✓ Sports Medicine,
- ✓ Metabolic Measurements,
- ✓ Human Performance laboratories,
- ✓ Physiological Human/Animal Research
- ✓ Flow Instrumentation
- ✓ Lung Simulation



HANS RUDOLPH, inc.
8325 Cole Parkway Shawnee, KS
66227, USA

PH: +1 (913) 422-7788
Canada: (800) 456-6695
FAX: +1(913) 422-3337
hri@rudolphkc.com
ISO 9001 - ISO 13485
www.rudolphkc.com


HANS RUDOLPH, inc.
International Marketing offices
PLM HealthCare Ltd
The Mills, Canal Street
Derby DE1 2RJ – UK



hri@plmhealthcare.com
PH: +33 (0) 612 181 501
FAX: +33 (0) 344 726 019

PLM Healthcare Ltd - V2™ FI APR 13 REV-1_9

Photos are not contractual. Specifications and features to change without prior notice.
© All rights reserved - Printed in the U.K.

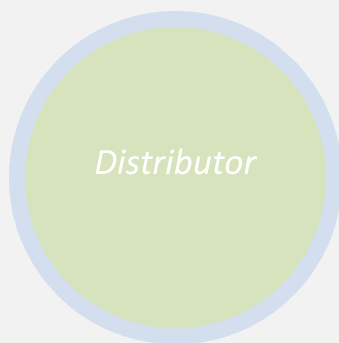

HANS RUDOLPH, inc.